

## AP5 pH/temp

## AP10 pH/mV/temp Features

Denver Instrument Portable Meters go where you go! Lightweight, ergonomically designed, waterproof meters with microprocessor accuracy

Four models in Denver's Portable pH meter lineup share a rugged, ergonomic design, with a sealed waterproof case that's buoyant and fully immersible in water up to 1 meter deep.

All models feature BNC and pin connectors, plus a waterproof Twist-Lock connector for special Denver Instrument electrodes.

Operation is exceptionally easy and the case fits comfortably in either hand, with soft rubber surfaces for a firm grip.

All measure pH with superb accuracy, with auto or manual temperature correction. And each one is ready to be toted most anywhere you go, providing up to 240 hours of use for Basic and about 40 hours for Advanced models on a set of four AA batteries. A selectable auto-off function saves batteries.




Two basic models offer accurate pH plus simple operation



- AP-5 for pH only
- AP-10 for pH, mV and relative mV readings

To operate, just press "Standardize" and the meter automatically recognizes any of 22 different buffers and checks electrode quality. You can calibrate with up to three buffers. Simple prompts and messages make accuracy easy.

Please contact Denver Instrument for more information:  
**1-800-321-1135 or 303-431-7255**  
or visit our website at [www.denverinstrument.com](http://www.denverinstrument.com)

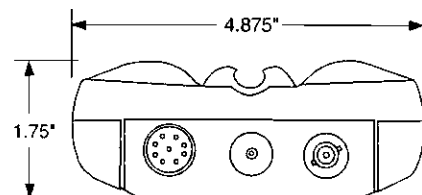

**Basic Portable pH Meter**

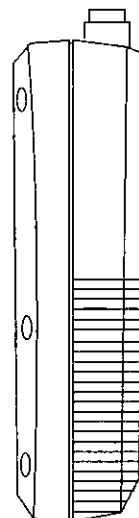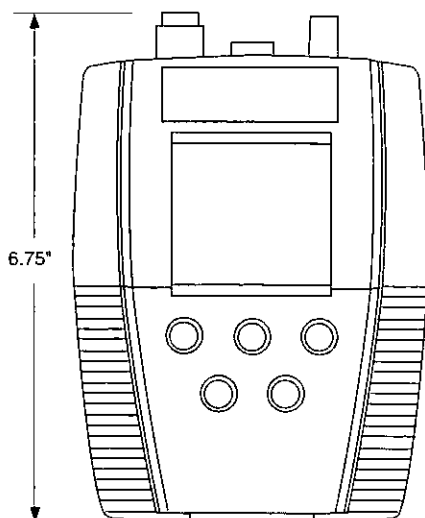

## AP5/AP10 Specifications

<b>pH</b>	Range	-1.99 to +19.99
	Resolution	0.01
	Accuracy	±0.01
<b>mV</b>	Range	-999 to +999 (AP10 only)
	Resolution	1
	Accuracy	±1mV
<b>Temperature</b>	Range	-5 to +105.0°C
	Resolution	0.1°C
	Accuracy	±0.4°C
<b>Standardization</b>	Up to 3 pH buffers manually or auto-recognize buffers from 4 sets (2,4,7,10,12: 1,3,6,8,10,13: 1.68, 4.0,6.86, 9.18, 12.46: 1.09, 3.06, 4.65, 6.79, 9.23, 12.75). 1 relative mV offset (AP10 only).	
<b>Dimensions</b>	4-7/8W x 6-5/8L x 1-3/4"H (12.5 x 17 x 4.5cm)	
<b>Power Requirements</b>	4 AA Batteries (included)	

**Meter with kit includes:** Meter, plastic-body pH/ATC "3-in-1" electrode, soft case, 4 AA batteries and operation manual.

**The Meter only includes:** Meter, 4 AA batteries, soft case and operation manual.

Please contact Denver Instrument for more information:  
**1-800-321-1135 or 303-431-7255**  
 or visit our website at [www.denverinstrument.com](http://www.denverinstrument.com)